Gebrauchsanweisung Operating Instructions Page mode d'emploi Istruzione per l'uso Instrucción

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Grasp the tickSAFE® Tickgripper in the palm of the hand with the thumb on the push button.

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
Press the push button completely down to fully open the gripper.

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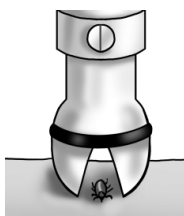Carefully position the opened gripper vertically to the skin surface so that the tick is positioned in the center of the two gripper tips.

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
Slowly release the push button. The gripper should grasp the tick and begin to rotate.

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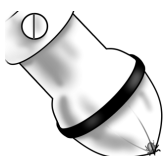The tick is removed and captured in the gripper. Should the removal procedure be unsuccessful, repeat Steps 1 to 4. Repeated attempts may be required for very small ticks!

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
Amtsgericht Nürnberg
HRB 25837, DE 268736024

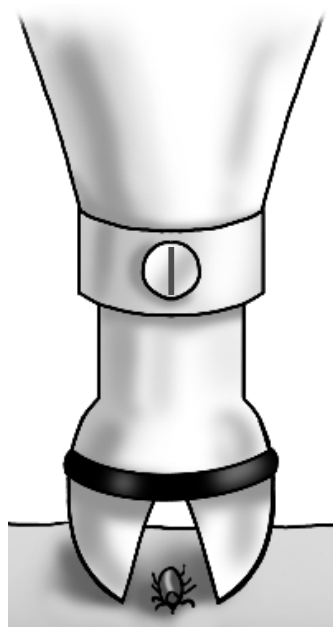Congratulations on the purchase of your new tickSAFE® Tickgripper! It will help you keep the risk of infection as low as possible when removing ticks. The unique patented technology utilized in the Tickgripper is proven effective. Please observe the following operating instructions:

1. Only use the Tickgripper to remove ticks! Do not misuse for other purposes!
2. As the Tickgripper contains small parts which a child may swallow, keep such safely stored away from small children.
3. Handle the Tickgripper carefully. Avoid flicking or playing with the push button. The internal mechanical parts can be damaged in such fashion. The Tickgripper is not a toy and thus unsuitable for small children.
4. Store the Tickgripper in such fashion that the gripper headpiece is not squeezed or compressed. Deformation can impair the precise closure of the gripper which is essential for the removal of small ticks. Replace the protective cap on the Tickgripper after use.
5. The screws on the gripper headpiece ring are pre-adjusted. Should with time a screw loosen, it can be carefully tightened. CAUTION: Overtightening the screw can cause the screw anchor to strip and/or hinder the proper rotation of the gripper.
6. Expose the Tickgripper to direct sunlight as little as possible. UV-rays can reduce the elasticity of the gripper headpiece material.
7. Large ticks can rupture during removal and contaminate the Tickgripper with blood. In such event the contaminated components can be cleaned with disinfectants. Frequent use of disinfectants can however reduce the elasticity of the gripper headpiece.

If you adhere to this advice, the Tickgripper will become an indispensable aid for tick removal. Our General Business Terms and Conditions can be found at www.ticksafe.com.

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